



ISBT 128 Blood Groups [ABO and RhD]

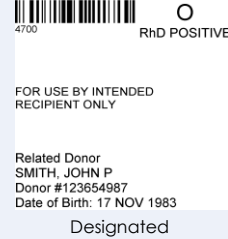

Questions:

		Yes	No	N/A	Comments
#1	Does the ABO/RhD appear in the upper right quadrant?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
#2	If autologous, directed, or designated donation, is the ABO/RhD printed BESIDE the corresponding bar code?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
#3	If the donation type is non-specified, is the ABO/RhD printed BENEATH the corresponding bar code?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
#4	Is the ISBT 128 ABO/RhD code printed left justified, immediately below, but not touching the linear bar code?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
#5	Does the ISBT 128 ABO/RhD code correspond to the ABO/RhD printed on the label?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
	Consult Table RT005 found in ST-001 .				
#6	Is the correct type of collection encoded in the ISBT 128 ABO/RhD code?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
	Consult Table RT005 found in ST-001 .				
#7	Is the type of collection encoded in the ISBT 128 ABO/RhD code non-conflicting with the type of collection encoded in the Product Code?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
	Consult Tables RT005 and RT008 found in ST-001 .				

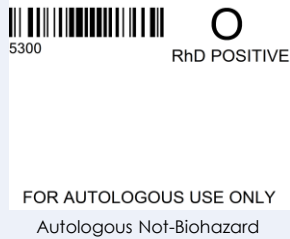

ST-001 ISBT 128 Standard Technical Specification
<https://www.iccbba.org/tech-library/iccbba-documents/technical-specification>

Audit Findings Summary:

Tally of responses

YES	NO	N/A
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Audit Tool for Cellular Therapy Version 1.0.1

© 2017 ICCBBA All Rights Reserved

Additional Auditor Questions:

a)

b)